# E>SWIN

FRANCE

**E-Flash 4c**

Français – English – Deutsch - Español

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## TABLE OF CONTENTS

<b>1. Prior to starting</b>	40
1.1 Removing hair at home using pulsed light technology	40
1.2 Overview of the <b>E-Flash 4C</b> device	41
1.3 Contents	42
1.4 Replaceable mobile head	42
<b>2. Intended use of the <b>E-Flash 4C</b> device</b>	42
<b>3. Safety with your epilator <b>E-Flash 4C</b> device</b>	43
<b>4. In which cases the <b>E-Flash 4C</b> device should not be used?</b>	46
4.1 Contraindications	46
4.2 Possible side effects with <b>E-Flash 4C</b>	47
4.3 Risks of inefficiency	50
<b>5. Precautions - How to use your <b>E-Flash 4C</b> device safely?</b>	51
5.1 Complying with safety and conditions of use	51
5.2 Using the device in a suitable environment	52
<b>6. Long-lasting hair removal method with <b>E-Flash 4C</b> device</b>	55
6.1 How does light manages to remove hair?	56
6.2 What is the impact of the hair growth cycle on photo-epilation?	56
6.3 Optimising results by planning hair removal sessions with <b>E-Flash 4C</b>	57
<b>7. Getting started with your hair removal device <b>E-Flash 4C</b></b>	58
7.1 Setting up your hair removal device <b>E-Flash 4C</b>	58
7.2 Removing your hair with the hair removal device <b>E-Flash 4C</b>	58
7.3 What should we expect with the hair removal device <b>E-Flash 4C</b> ?	61
7.4 After treatment with your hair removal device <b>E-Flash 4C</b>	64
<b>8. Maintaining your <b>E-Flash 4C</b> device</b>	62
8.1 Maintenance and cleaning of the <b>E-Flash 4C</b> device	62
8.2 Changing the replaceable mobile head	63
<b>9. Troubleshooting</b>	65
<b>10. Customer Service</b>	66
<b>11. Frequently Asked Questions</b>	66
<b>12. Specifications</b>	71
<b>13. Packaging - Labelling</b>	72
<b>14. Warranty</b>	72

## 1. Prior to starting

### 1.1 Prior to starting

The **E>Flash** range which includes the device you have purchased has been developed with all the know-how of French company **E>Swin**, designer and manufacturer of high-tech medical devices.

Four years have been required to develop a flash containing amazing features: having the same strength as the devices of competing big brand, its effectiveness in terms of hair removal has been greatly increased.

Clinical studies conducted by an independent laboratory have confronted the **E>Flash** to big brand products. The results are final and the difference is visible from 15 days after the first session.

In parallel, we have put all our skills in developing manufacturing procedures for **E>Flash**, for the purpose of ensuring full traceability and consistency of integration of all electronic components.

Finally, we hope that your **E>Flash** will bring you the full satisfaction that you are entitled to expect from such a product. You can discover all the fun of a modern hair removal with surprising results as from the first sessions.

Yours sincerely,

*Yves Vincent*

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## 1.2 Overview of the E-Flash 4C device

The hair removal device E-Flash is a long-lasting photodepilation device that can be used in the comfort of your home.

Your device has three power levels, a skin colour detector and a replaceable mobile head located on the tip of the device, as follows:

:

## 1.3 Contents

- 1 hair removal device **E-Flash 4C** (including replaceable mobile head enabling 5000 flashes)
- mains adapter
- 1 instruction manual

## 1.4 Replaceable mobile head

The mobile head of the hair removal device **E-Flash 4C** contains a lamp with a capacity of 5000 flashes (whatever the level of intensity). Depending on your usage and the number of flashed areas, you may need to change the replaceable mobile head.

## 2. Intended use of the **E-Flash 4C** device

It may be that the use of the hair removal device **E-Flash 4C** at home is a new experience for many people in terms of hair removal by light. The hair removal device **E-Flash 4C** is easy to handle and hair removal sessions are done quickly.

Removing hair with **E-Flash 4C** device enables people to obtain lasting and visible results as early as 15 days after the first session.

The hair removal device **E-Flash 4C** is a hair removal device for removing unwanted hair from the body, including legs, arms, underarms and bikini lines.

The **E-Flash 4C** device is for adult men and women, as well as minors under the supervision of an adult.

### 3. Safety with your epilator E-Flash 4C

Safety is a priority when using the hair removal device E-Flash 4C; the device's technology offers maximum safety and minimum risk.

#### REMINDERS OF THE SAFETY INSTRUCTIONS RELATED TO LUMINOUS POWER

##### Power level Intensity

The power level determines the power of pulsed light administered on the skin during treatment, from lowest to highest. The corresponding blue power level indicator light indicates the level of power selected. By increasing the power, you get better results, but you also increase the risk of side effects (see § 4.2 "Possible Side Effects with E-Flash 4C").

Always use the minimum power the first time you use your hair removal device E-Flash 4C.

If during your first hair removal session or after, you feel little or no discomfort, you could try, for your next session, to increase the power level by a notch, which will be indicated by the corresponding blue indicator light and so on.

##### Do not increase the power:

- If you feel more discomfort than usual after or while using the device (see § 4.2 "Possible side effects with E-Flash 4C")
- If you do not feel any discomfort during the session. It is possible that you do not feel anything during treatment, but you might feel some discomfort afterwards.

When your device is turned on, the intensity level is automatically set by default to the lowest intensity. The blue indicator light below the power indicator will be lit. To adjust the light to the desired intensity, either higher or lower, press the power selector switch until the indicator light that corresponds to the intensity turns on.

### The skin colour detector

The hair removal device **E-Flash 4C** gets rid of unwanted hair by targeting the melanin in the hair. The amount of pigment in a person's skins determines its colour, but also the level of risk to which it is exposed when using a method of hair removal by light, regardless of the method. This is why hair removal based on pulsed light is not convenient for naturally dark skins and carries a risk of side effects, including heating of the skin, burns, phlyctenes, and skin discoloration [hyper or hypo-pigmentation].

The hair removal device **E-Flash 4C** was designed with a skin detector device that enables it to check the colour of the skin to be treated before each flash. This sensor verifies that the energy delivered is only administered to suitable skin types. This device will not allow you to treat areas where the skin is too dark or too tanned.

Also, the optical window through which the pulse of energy is issued, and which represents a surface area of 4 cm<sup>2</sup>, is embedded in the mobile head. These elements will provide optimum protection by preventing direct skin contact with the optical window.

**Do not use the hair removal device **E>Flash 4C** on dark skin, it will not work.**

- Do not use the hair removal device **E>Flash 4C** on tanned skin or after recent sun exposure.
- Avoid exposure to sunlight during the 4 weeks that precedes depilation with the **E>Flash 4C** device and for 2 weeks after flashing. "Exposure to sunlight" means an unprotected exposure to indirect rays during a continuous hour. You must protect your skin with sunscreen having a sun protection factor of 50 or more, throughout the entire treatment period and at least 2 weeks after the last hair removal session.
- Do not expose treated areas to sunlight.

#### **Eye protection device**

- Flashing on skin can cause serious eye damage.
- **E>Flash 4C** eye protection device prevents all flashes from being issued if device is pointed towards an empty space. **E>Flash 4C**, in fact, can only work when it is in contact with and pressed onto the skin.

#### **4. In which cases the E-Flash 4C device should not be used?**

Do not use the E-Flash 4C device in any of the following situations :

- Pregnant women
- People with diabetes.
- Dark skins, tanned skins (of phototype V and VI, or even IV).
- Skin that has been highly exposed to sunlight or UV machine during the 28 days that precedes treatment (see Chapter 3).
- Tattoo or permanent makeup on the area being treated.
- Wounded, burned, or infected skin.
- Skin lesions with unusual appearance (medical opinion should be sought).
- Moles (if they are located on the area to be treated, they must first be hidden), brown or black lesions, pregnancy marks, freckles, birthmarks and pigmented nevus.
- Eczema, psoriasis, lesions, open wounds or active infections on the area to be treated. Wait until the affected area is healed before using the E-Flash 4C device.
- Treatment for photosensitivity (medical advice should be sought).
- When there is a history of skin disease.
- Abnormal skin caused by diabetes, for example, or by other systemic or metabolic diseases.
- All other cases that your doctor may consider risky for treatment.
- Do not use the hair removal device E-Flash 4C on the nipples and genitals (male or female).
- Do not use the hair removal device E-Flash 4C other than for hair removal.
- Do not use on animals.

In case of doubt about the personal use of the hair removal device E-Flash 4C, please consult your doctor or dermatologist.

#### 4.2 Possible side effects with E-Flash 4C

When used properly, the risk of adverse effects, such skin rashes or burns, is extremely reduced but still exist. If these effects occur, you must stop using the device immediately.

A simple solution is to test the device on a small area of the skin before the full session. The absence of skin reaction will enable a perfectly safe use of the device.

When your hair removal device E-Flash 4C is used as directed, side effects and complications are rare. However, all cosmetic treatments, even those designed for home use, involve certain risks. This is why it is important to understand and accept these risks and possible complications when using a pulsed light hair removal system for home use.

Avoid overlapping flashes

**In case of overlapping or repeated administration of pulses on the same area:**

Side Effects	Demonstration of side effects
Discomfort / unpleasantness on skin	Light
Slight redness	Light
Skin hypersensitivity	Light
Burns or wounds	Rare
Scarring	Rare
Change in pigmentation	Rare
Excessive skin redness and swelling	Rare
Infection	Insignificant
Bruising	Insignificant

**• Feeling of a mild discomfort on skin**

Although hair removal by pulsed light at home is generally well tolerated, most people who use it feel a slight discomfort during use often described as a tingling sensation on the treated area. In general, this feeling lasts the length of the application, maybe for a few minutes afterwards. Any discomfort beyond this unpleasantness is abnormal and means that you should either stop treatment because you cannot tolerate it, or that the power level you have selected is too high.

- **Redness**

Your skin may become red immediately after using the hair removal device **E-Flash 4C** or within 24 hours. The redness usually disappears within 24 hours. Consult your doctor if it persists beyond 2 to 3 days.

- **Excessive redness and swelling**

In some rare cases, the treated skin may become very red and swollen. This reaction is more common on sensitive body parts. The redness and swelling should disappear within 2-7 days during which you will need to apply ice regularly. A gentle cleaning does not pose any problem but you should avoid sun exposure.

- **Skin hypersensitivity**

Skin is sensitive on the treated area and could become dry or flaky.

- **Change in pigmentation**

The **E-Flash 4C** device targets the root of the hair, especially the melanin in the hair bulb. However, there is a risk of temporary hyperpigmentation (increased pigmentation or brown discoloration) or hypopigmentation (whitening) of the surrounding skin area. The risk of changes in pigmentation is higher among people with dark skin.

Generally, the change in pigmentation and discoloration are temporary. Cases of permanent hyperpigmentation or hypopigmentation are rare.

- **Scarring**

Although these cases are very rare, permanent scarring can occur. Generally, scars take the form of a superficial whitish lesion (hypotrophic). Scars can be large and red (hypertrophic) or large and spread beyond the injured area (keloid). More advanced aesthetic treatments may be required to improve its appearance.

- **Burns or wounds**

In very rare cases, the skin can be exposed to injuries or burns following the application. Several weeks may be required for a full recovery and, in exceptional cases a scar may remain.

- **Infection**

Skin infection is extremely rare, but can occur after a burn or skin damage caused by the hair removal device **E>Flash 4C**.

- **Bruising**

In very rare cases, the use of the hair removal device **E>Flash 4C** can cause a blue purple bruising that can last between 5-10 days. As the blue colour disappears, the skin may take on a rust colour (hyperpigmentation) that is potentially permanent.

Stop treatment immediately if skin shows signs of burns or blisters.

For more information or personal advice, please visit [www.e-swin.com](http://www.e-swin.com) or your area's approved distributor's website..

#### 4.3 Risks of inefficiency

The use of **E>Flash 4C** is not suitable for everyone. The **E>Flash 4C** device is not efficient for the treatment of hair of the following colours: white, grey, red and light auburn..

- The **E>Flash 4C** device is not effective on very light, almost white hair.

## 5. Precautions - How to use your E-Flash 4C device safely?

E-Flash 4C is a long-lasting hair removal device using pulsed light professional technology. As such, it is not a trivial device and therefore complying with the precautions for usage and restrictions of use is essential.

### 5.1 Complying with safety and conditions of use

We particularly recommend:

- Selecting the intensity level of power with caution
- Avoiding side effects.

To do this, comply with the instructions for use and the contraindications described in Chapters 3 and 4.

## 5.2 Using the device in a suitable environment

This appliance is intended for household use only.

As with any electrical equipment, certain precautions are necessary to ensure your safety. Before using the **E>Flash 4C**, please read the following information carefully:

- Voltage: 100-240 V ~.
- Nominal frequency: 50-60 hz.
- Maximum power consumption: 540 VA.
- Make sure that the voltage corresponds to that of your electrical installation. Improper connection may cause permanent damage not covered by warranty.
- Make sure that the power outlet to which your **E>Flash 4C** is connected is undamaged, in perfect working condition and that it does not interfere with the proper use of your **E>Flash 4C** device (the power outlet must be easily accessible). Do not connect the hair removal device **E>Flash 4C** if the adapter cord or plug is damaged.
- Do not use any extension cord or cable or adapter other than those provided by E-Swin as this could render the equipment non-compliant on an electromagnetic point of view.
- Do not stack your device on other electronic devices.
- Never leave your hair removal device **E>Flash 4C** unattended while switched on.
- Always unplug the adapter of the **E>Flash 4C** device from the mains supply immediately after use, without ever pulling on the cord of the **E>Flash 4C**.
- Unplug your **E>Flash 4C** when cleaning. Do not use or stop using the device if you notice some anomaly, or if it shows signs of defect or make a strange noise, or has a strange odour. More particularly, cease using it immediately if the replaceable mobile head, the anti-UV yellow protection filter or skin colour detector show signs of a crack or stains.
- Do not use the **E>Flash 4C** following a fall (non-visible damage to the device may be detrimental to your safety).

- Please contact the customer service of your hair removal device **E-Flash 4C** if your device is broken, damaged or requires repair.
- Do not use accessories other than those recommended by **E-Swin**.
- E-Swin disclaims any liability for any damage resulting from improper use of the device or use other than that described in this manual.
- Never use the **E-Flash 4C** near an aerosol, a heat source, or near other electronic devices.
- The wireless communication equipment must be kept away (about 3 meters) because they may interfere with proper functioning of the device.
- The **E-Flash 4C** must be placed in an environment relatively unpolluted by electromagnetic interference (computer, electrical appliances...).
- Do not use accessories or components with your **E-Flash 4C** device that have not been recommended by **E-Swin**. **E-Swin** disclaims any responsibility for the use of replaceable mobile heads and accessories from another manufacturer.
- When **E-Flash 4C** is running, it should not come into contact with hair, face, clothing or other objects, so as to prevent all risks of injury, blockage or damage.
- Never use the **E-Flash 4C** with wet hands.
- **E-Flash 4C** must not be used in a place where there is a possibility for water to splash the device.
- Do not store the hair removal device **E-Flash 4C**, the replaceable mobile head or the adapter in a place where it can fall or find itself in a bathtub or sink, or any other area containing liquid. Do not place in water or any other liquid. This will cause a serious electrical shock.
- Do not use the hair removal device **E-Flash 4C** in your bath or shower.
- Unplug the **E-Flash 4C** device immediately if it gets dropped in water.
- Never look directly at the light coming from the replaceable mobile head.
- Never point the application surface of the hair removal device **E-Flash 4C** on an open space to attempt to make a flash. Make sure that you always point the treatment surface to the skin and that it is fully in contact with the skin before releasing a flash.

- Keep the device out of reach of children. Do not allow children to use the hair removal device **E>Flash 4C**.
- This device must not be used by people (including children) with reduced physical, sensory or mental capabilities, or by people deprived of experience or knowledge, unless they have been supervised or have been given instructions by a responsible guardian concerning the use of the device. Children should be monitored so that they do not play with the device.
- Never use flammable liquids such as alcohol (including perfume, deodorant and other composition containing alcohol) or acetone to clean the skin before using the hair removal device **E>Flash 4C**. Using your **E>Flash 4C** device can cause a temporary change in pigmentation (see § 4.2 "Possible Side Effects with **E>Flash 4C**."

**Keep the hair removal device **E>Flash 4C** away from any water sources.**

When using the device in a bathroom, be sure to disconnect it after using it. Indeed, being near a source of water can be dangerous even if the power is switched off. To provide additional protection, installing a residual current device (RCD) of a maximum rating of 30mA inside the electrical circuit supplying the bathroom, is recommended. Ask your installer for advice.

Warning! The polyethylene bags containing the product or its packaging can be dangerous. Keep these bags away from babies and children. Do not use in cribs, cribs, strollers or playpens. The thin film may cling to nose and mouth and prevent breathing. A bag is not a toy.

**Do not attempt to open, modify, or repair the hair removal device E-Flash 4C on your own. Only repair centres approved by E-SWIN are authorized to make repairs.**

Trying to open the hair removal device E-Flash 4C may damage the device and void your warranty.

### **End-of-life electrical and electronic equipment**

In the interest of all and to actively participate in the collective effort of environmental protection: Do not discard your product with household waste. Use return and collection systems that are available in your area. Some materials may be recycled or recovered.

Note: The E-Swin showroom located in Paris (France) in the 17th arrondissement (34 rue Brunel) collects used mobile heads.

## **6. Long-lasting hair removal method with E-Flash 4C device**

The hair removal device E-Flash 4C is a long-lasting personal photo-epilation device. The process of hair removal by pulsed light is well known and understood. It has been clinically tested around the world for over 15 years and has proven its effectiveness and safety to offer gradual long-term hair removal.

## 6.1 How does light manage to remove hair?


**E-Flash 4C** uses pulse light technology that enables slower hair regrowth and reduces hairiness gradually as and when hair removal sessions happen.

Photo-epilation is based on the theory of selective photothermolysis using optical energy to stop hair growth. In order to achieve such thermal effect, the hair root must absorb photo energy in a selective way and convert it into heat. This selection occurs when the optical energy issued is mainly absorbed by the pigment of the hair bulb, responsible for the absorption of light, which generates the heat that will stop hair growth. Also, the more melanin the hair contains (i.e. the darker the hair is), the more light it can absorb and thus, photo-epilation becomes effective.

## 6.2 What is the impact of the hair growth cycle on photo-epilation?

Every hair in our body goes through three phases of growth: the anagen phase, the catagen phase and the telogen phase. These three phases play an important role in the photo-epilation process. The anagen phase is the growth phase of the hair, while the catagen and telogen phases are both resting phases. The time it takes for a hair to get through its entire hair cycle varies from person to person, and depends on the region of the body where it is located. However, the average is between 18 and 24 months.

The majority of hair follicles found themselves in one of these two resting phases at one point. Therefore, the hair removal device **E-Flash 4C** will have no effect on the hairs when they are in these resting phases. However, it will act on hairs that are in the anagen phase
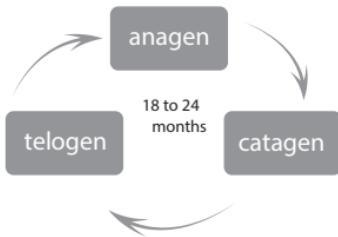
## 6.3 Optimising results by planning hair removal sessions with **E-Flash 4C**

Care program with **E-Flash 4C** throughout a full hair growth cycle:

- The first 4 sessions with the hair removal device **E-Flash 4C** must be spaced at a two weeks interval. (Over a period of about 2 months).
- The fifth session and subsequent sessions must be spaced at one-month interval.
- Over time, you will be able to gradually increase the intervals between the sessions according to the results obtained.

Given the hair growth cycle, several sessions with your hair removal device **E-Flash 4C** will be needed to achieve sustainable results. The effectiveness of the device varies from one person to another depending on the area to be treated, the type of hair and biological factors that can determine the hair growth cycle area. How you use your product is also a determining factor.

Some people will get results faster or slower than the average number of sessions that we recommend. One cannot expect lasting results in just one or two sessions.

You can plan your personal sessions at your own pace and still get satisfactory results. However, we recommend that you follow the schedule described above because it has shown excellent results.

Note: treating the same area more than once in two weeks does not give better results and may cause adverse effects.

Due to hormonal factors and physiological changes, resting follicles can become active. Maintenance sessions may therefore be necessary from time to time, even after you have reached your goal.

## **7. Getting started with your hair removal device E-Flash 4C**

### **7.1 Setting up your hair removal device E-Flash 4C**

1. Remove the device and all other components of the hair removal device E-Flash 4C from the box..
2. Make sure that the replaceable mobile head is properly inserted at the tip of the applicator.
3. Plug the adapter cord into the slot provided for this purpose at the rear of the device.
4. Connect the other end of the adapter into an electrical outlet.

Your hair removal device E-Flash 4C is ready to use.

### **7.2 Removing your hair with the hair removal device E-Flash 4C**

The representations mentioned in the following paragraphs are shown on page 5.

#### **Prior to each use:**

Make sure to check your device in its entirety, including the external optical parts of the E-Flash 4C .

## Step 1: Prepare your skin

- The skin should be clean, dry and free of any cosmetic product.
- Make sure that no cream, moisturizer, oil, or any other substance has been previously applied to the skin (including during the 12 hours prior to the session).
- Always shave the area to be treated with running water, on the day of the session or one hour before the session for the most sensitive skin (do not use any methods other than shaving to remove the hair, as it would make the treatment with pulsed light inefficient).
- Dry the skin (and do not apply any substance on the skin).

A first test session on a small area of the skin the day before the full session will allow you to make sure that your skin supports the treatment. The absence of skin reaction will enable a perfectly safe use of the device.

## Step 2: Prepare your E-Flash 4C

- Once the E-Flash 4C is connected (see section 7.1), hold down the “on/off” button located at the rear of the device **C**. for about 2 seconds. Subsequently, the three blue lights start flashing and the fan starts working (**i**, **j** et **k**).
- When the device is ready to be used, the first blue indicator light remains lit (**i**).
- Choose the power level intensity (1, 2 or 3) of the flash, from the bottom to the top of the device:
  - The 1st blue indicator light corresponds to power level 1 (**i**)
  - The 2nd blue indicator light corresponds to power level 2 (**j**)
  - The 3rd blue indicator light corresponds to power level 3 (**k**)
- To go to power level 2 (represented by the light **i**), press once on the button located on the rear of device **C**.
- To go to power level 2 (represented by the light **k**), press again on button located on the rear of device **C**.
- To reduce the power of the device, press the button located on the rear of the device until you reach the desired power level.

### **Step 3: Flashing**

1. Grab the device, making sure that your hand does not obstruct the air vents (found under the device and its rear tip).
2. Press the mobile head of the device to the skin; making sure that its entire surface is fully in contact with the skin.
3. Hold this position and trigger the flash by pressing simultaneously on the two buttons found on either side of your device (right and left , **a** and **b**).

**Please refer to Chapter 3 of this manual to determine the power level intensity.**

The more tanned or naturally pigmented the skin is, the more you must reduce the power level to avoid any overheating or burning phenomenon.

If you do not see any flash of light and the red light turns on, it means that the skin is too dark for safe application. Change the treatment area or check with the help service of the hair removal device **E-Flash 4C** (see details 3, "Skin colour detector").

Your **E-Flash 4C** device will charge immediately between each flash, the delay between the onset of each flash is up to 2.9 seconds. During this period, the first blue indicator light indicating the power level will turn on, and then successively, the second will turn on and finally the third one will as well. When the device is ready to flash, the blue indicator light(s), depending on the power level you have chosen, will stay fixed.

4. Move the mobile head from the application surface of treated area.
5. Move the device to the next area and start flashing. You must administer the pluses in a row, from one end to the other, and by proceeding methodically through **successive rows**. The mark that the application pressure leaves on your skin should help you to align the device correctly before triggering a new flash. This method helps you to better cover the area in question while preventing you to overlap the flashes or to issue them too close to each other on the same area, or on the other hand, to leave too much space in between the areas.

**Do not treat the same area more than once per session**, otherwise the risk of side effects will increase (see § 4.2 "Side Effects with **E-Flash 4C** »).

Please note: If the **E-Flash 4C** device stays idle for a prolonged period when in used, it will automatically turn off, you will once again need to hold down the "on/off" button, located at rear of the device, for about 2 seconds". Also, when the power is turned on, if you hold down the "on/off" button for too long (about 6 seconds), the **E-Flash 4C** Flash will turn off again.

## 7.3 What should we expect with the hair removal device **E-Flash 4C** ?

Here are normal things that you will see or hear during a session:

- **Fan noise**

The cooling fan of the hair removal device **E-Flash 4C**, makes a noise resembling that of a hair dryer, this is normal. **E-Flash 4C**

- **Soft sound at each time it issues a flash**

Upon activation of a new flash of the hair removal device **E-Flash 4C**, a soft sound will accompany the luminous flash.

- **Moderate pressure from the replaceable mobile head**

This is necessary and useful to correctly position the device on the area to be treated; this is exclusive to the safety device of hair removal device **E-Flash 4C**.

- **Tingling and sensation of heat**

It is normal to feel a moderate sensation of heat and tingling caused by the light energy at each flash. Remember, it is important to use the minimum power levels for the first sessions. The heat sensation can last up to 1 hour.

Other possible side effects are also possible, please refer to § 4.2 "Possible Side Effects with **E-Flash 4C**".

## 7.4 After treatment with your hair removal device **E-Flash 4C**

- At the end of your session with your hair removal device **E-Flash 4C**, switch off your device by holding down the "on/off" button for 3 seconds. Disconnect the adapter from the power outlet.
- We recommend that you clean your hair removal device **E-Flash 4C** after each application session (see "Care and of the **E-Flash 4C** device").
- After cleaning your device, we recommend storing it in its box and keep it away from any water source.

## 8. Maintaining your **E-Flash 4C** device

### 8.1 Maintenance and cleaning of the **E-Flash 4C** device

We recommend that you clean your hair removal device **E-Flash 4C** after each session. Unplug the device before cleaning.

- **Plastic parts**

All plastic parts of the device must be cleaned with a clean cloth, dampened with water (without any additives). Be careful not to drop any liquid into the air vents located at

the rear and at the tip of the device: they must be kept clean without elements likely to obstruct them.

- **Power supply**

Check that the power supply and the cord are in perfect condition (must not be pinched or crushed, no bare wires). Clean with a cloth dampened with water (without any additives).

- **Mobile head and application surface**

The external part (especially the anti-UV yellow protection filter) can be cleaned with a dry cloth or with a clean and dry cotton swab.

Never immerse your **E-Flash 4C** device or one of its components in water.

## 8.2      Changing the replaceable mobile head

**E-Flash 4C** is equipped with a replaceable mobile head containing a lamp that issues 5000 flashes.

The intensity of the flash is defined by adjusting the power level of the device. There is no loss of energy in the life of a replaceable mobile head. When there are only 500 flashes, then 250 flashes, then 100 flashes successively left on the mobile head, the red indicator light will warn you (see Chapter 9, "Solving problems").

When the 5000 flashes of the replaceable mobile head have been issued, the red indicator light that indicates that the mobile head needs to be replaced will turn on: this has been 100% used. This light warns you that no other flash will be issued. To continue with your hair removal process, you should replace the replaceable mobile head.

To change the replaceable mobile head:

1. Hold down the "on/off" button of your hair removal device **E-Flash 4C** for about 3 seconds.
2. Unplug the adapter from the power outlet.
3. Remove the used mobile head by pressing simultaneously on the two locking tabs located on either side of the mobile head (1) and then release them.

A small "click" means that the head can now be removed. To dispose of the used head, use the return and collection systems that may be available in your area or leave it at our **E-Swin** showroom at 34 Brunel road in Paris (75017).

4. Unpack a new mobile head.

Note: the mobile head must be replaced prematurely if large spots appear on its lamp or inside or if its mobile head, its application surface or anti-UV yellow protection filter are damaged.

## 9. Troubleshootings

Problem	Solution
My device does not start (it does not turn on and the indicator lights are not working)	<p>1. Make sure the adapter is well connected to your hair removal device <b>E&gt;Flash 4C</b> and connected to power outlet.</p> <p>2. When you press the "on / off" button of the device, make sure you are not holding it down for more than 6 seconds.</p>
The red indicator light flashes 2 times in a row (every 2 seconds) and the blue indicator lights are working.	There are less than 500 flashes left on the cartridge included in the mobile head
The red indicator light flashes very quickly (every 0.4 seconds) and the blue indicator lights are working	There are less than 250 flashes left on the cartridge included in the mobile head
The red indicator light stays on and the blue indicator lights are working	There are less than 100 flashes on the cartridge included in the mobile head
There is no flash when I press the trigger buttons (the red light stays on but the device is not issuing any flash)	There are no more flashes on the cartridge included in the mobile head
Il n'y a pas de flash lorsque j'appuie sur les boutons déclencheurs (le voyant rouge reste allumé mais l'appareil refuse de flasher).	<p>Make sure that the mobile head is positioned evenly and is firmly in contact with the skin. For your safety, the trigger will be able to issue a flash only if the moving head is pressed firmly against the skin.</p> <p>Make sure that the mobile head is securely attached to your device.</p> <p>Make sure that the skin to be treated is not too dark. Try using the device on another area.</p> <p>Make sure that the mobile head has not run out of number of flashes.</p> <p>If unsuccessful, please contact our customer service department.</p>

## **10. Customer Service**

For more information on **E-Flash 4C** products, please visit **E-Swin**'s website on [www.e-swin.com](http://www.e-swin.com) site.

If your device is broken, damaged, needs repair or if you need assistance, please contact **E-Swin**'s customer service department.

## **11. Frequently Asked Questions**

### **1. Is the hair removal device **E-Flash 4C** really efficient?**

Yes. Clinical studies conducted by an independent laboratory have confronted the **E-Flash 4C** to big brand products. The results are final and the difference is visible as early as 15 days after the first session.

Also, **E-Swin**, which is the creator of the hair removal device **E-Flash 4C** specialises in pulsed light technology for various applications that are dedicated to professionals as well as the general public.

### **2. On which body parts can I use the hair removal device **E-Flash 4C** ?**

The **E-Flash 4C** device has been designed for the removal of body hair. The areas most commonly treated by the hair removal device **E-Flash 4C** are legs, arms, armpits and bikini line. The **E-Flash 4C** device cannot be used on the face.

### **3. How long lasts a hair removal session with hair removal device **E-Flash 4C**?**

The duration of treatment depends on the area being treated. Both lower legs may take on average 20 minutes and both armpits 2 minutes. The hair removal device **E-Flash 4C**

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#### **4. Is the hair removal device **E>Flash 4C** safe?**

Your safety is our priority. When we designed the hair removal device **E>Flash 4C**, it was tested and approved according to strict standards for home use. Like any electronic device in contact with the skin, the user must follow the instructions carefully and comply with the safety instructions.

#### **5. Are treatments with the hair removal device **E>Flash 4C** painful?**

When used properly, most users of the hair removal device **E>Flash 4C** reported a slight sensation of heat when the flash is issued. Users who have darker and thicker hair feel a little discomfort, but these sensations subside and disappear at the end of the session. For your convenience, the hair removal device **E>Flash 4C** offers 3 power levels depending on your sensitivity.

#### **6. How often should I use the hair removal device **E>Flash 4C** ?**

The first 4 sessions with hair removal device **E>Flash 4C** must be spaced at a two weeks interval. The following sessions must be spaced at a four weeks interval, and then the sessions will be carried out as the hair grows back until desired results are achieved.

#### **7. Is the hair removal device **E>Flash 4C** effective on white, grey or light auburn hair?**

The hair removal device **E>Flash 4C** is most effective on dark hair or those containing more melanin. Melanin, the pigment that gives colour to the skin and hair (and hair) absorbs optical energy. Black, brown, and dark hair responds best to treatment. If brown and light brown hair also responds to treatment, they require a larger number of sessions. White, grey, red, or light auburn hairs do not respond to treatment.

## **8. Can I use the hair removal device E-Flash ЧС on a dark or black skin?**

Do not use the device on naturally dark skin. The hair removal device E-Flash ЧС removes hair by targeting hair pigment. Pigments are also found in the surrounding dermal tissue in various quantities. The amount of pigment in the dermis of a person visible by his skin colour makes it possible to define the level of risk to which he is being exposed while using the hair removal device E-Flash ЧС. Treating dark skin with the hair removal device E-Flash ЧС may cause burns, blisters or a change in skin colour (hyper or hypo-pigmentation). An exclusive safety device - a skin colour detector - has been designed inside the hair removal device E-Flash ЧС to check the colour of the skin at the beginning of each session and occasionally during the current session. When the device detects a darker area, the sensor disables the hair removal device E-Flash ЧС which then cannot issue any flash.

## **9. When will I start seeing results with the hair removal device E-Flash ЧС ?**

As with any other method of hair removal by the light (laser or pulsed light), the results are not instant; and you might even have the impression that nothing has happened. The hair sometimes gives the impression of continuing to grow after a session, but after two weeks, most of these hairs will simply fall off. Moreover, the hair growth is done in 3 phases and only hairs found in the active phase are affected by the hair removal device E-Flash ЧС; which is why it is necessary to multiply the sessions to achieve the expected result.

## **10. Can a man use the hair removal device E-Flash ЧС ?**

The hair removal device E-Flash ЧС is suitable for men, but male hairs, especially on the torso, require a greater number of sessions to achieve the desired results.

## **11. Why does hair grow back on the area that I treated last week?**

It is quite common that the hair gives the impression of continuing to grow for 1-2 weeks after a hair removal session with the hair removal device **E-Flash 4C**. This process is called "ejection." After 2 weeks, however, you will find that these hairs are falling off or are releasing themselves from their follicle. That being said, we do not recommend that you pull on the hairs to release them. Let them fall naturally. In addition, some hairs will not have been affected by the hair removal device **E-Flash 4C** due to either misapplication or the fact that the hairs are in a resting phase. These hairs will be treated during subsequent sessions, hence the need to have several sessions in order to achieve the desired effect with the hair removal device **E-Flash 4C**.

## **12. I heard that some hairs grow back lighter and finer after photo-epilation?**

This phenomenon is well documented by professionals who use devices with light technology (Laser and pulsed light) for hair removal treatments. It is possible that some hairs grow back lighter and finer after depilation with the hair removal device **E-Flash 4C**. In general, these hairs are only a small part of the original hairiness. Treatment must be continued to achieve the desired result.

## **13. Is long-term use of the hair removal device **E-Flash 4C** dangerous for my skin?**

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**14. Can I use the hair removal device E>Flash 4C to remove hair from the chin or face?**

The hair removal device E>Flash 4C from E>Swin cannot be used for the treatment of facial hair.

**15. How to maintain the treated area after using the hair removal device?**

The area treated with the hair removal device E>Flash 4C can be cleaned and maintained with recognised and hypoallergenic products for skin care. Special precautions must be taken to avoid unprotected sun exposure. Sunscreens with high protection factor (over 50) and long clothing are a good protection against the sun.

**16. Should I pluck the hairs before, between and after the treatment?**

No, do not pluck the hairs before, between and after treatments. The hair will fall off naturally and gradually. This stage can last up to 2 weeks. If you wish, you can shave the treated areas between sessions (but do not pluck or pull hair).

## 12. Specifications

Impact point, specifications	1,33 cm x 3 cm (4 cm <sup>2</sup> )	
Speed	1 flash every 2.5 seconds	
Technology	Pulsed Light	
Light intensity	3,5 to 5 J/cm <sup>2</sup>	
Wavelength	470 - 1200 nm	
Charging time / energy source	Continuous operation	
Electrical Installation Requirements	100 – 240 VAC, 0,9 A	
Time required for the treatment of lower legs	20 minutes	
Handling and safety	A skin colour detector enables safe use	
Temperature	While in operation	10 to 25° C
	Storage	- 10 to + 50 ° C
Relative Humidity	While in operation	30 to 75 %
	Storage	10 to 90 %
Atmospheric pressure	While in operation	700 to 1060 hPa
	Storage	50 to 1060 hPa

## 13. Packaging - Labelling

Degree of protection against electric shock: BF Type device

Waste of Electrical and Electronic Equipment

CE Marking

Follow directions for use

Protection against moisture: Normal

This product must not be used in the presence of a flammable anaesthetic mixture with air or with oxygen or nitrous oxide.

## 14. Warranty

### **What does this warranty cover?**

**E-Jwin**, or the approved distributor from whom you purchased your **E-Flash 4C**, (hereinafter referred as the Device) will repair or replace the Device, if any hardware defect caused by **E-Jwin** is found in the device during the warranty period, whether it is a hardware defect on the body of the device or a malfunction.

### **Who is entitled to the services covered by this warranty certificate?**

Only you, the original purchaser of the Device, who have purchased the device from **E-Jwin** or its approved distributors, are entitled to benefit from the warranty specified on this certificate. To benefit from the warranty service, you must prove that you are the original purchaser by showing the original invoice of the device. This warranty is not

an international warranty. It will be honoured only in the country in which the approved distributor or the **E-Swin** entity from whom you bought the device is located.

### **What is the warranty period covered by this certificate?**

The warranty covers repair requests received (a) within one (1) year from the date of purchase of the Device or (b) ninety (90) days from the date on which the repaired or replaced Device has been handed to you, whichever is the latest. If a self-adhesive label, indicating a warranty expiry date that is further away, is attached to your Device, the warranty expiry date stated on the label is the one that applies, and the device will be repaired or replaced under this warranty until the date indicated.

If you have purchased the device in a European Union country, you get an extra year (in accordance with EU Directive 1999/44) to present a claim for defect or malfunction of your device, which exist but not apparent at the time of purchase, as determined by **E-Swin** or its approved distributor. For this warranty extension to be effective, you must contact the approved distributor from whom you have purchased the device or local repair service and notify the production defect without delay and in within a maximum period of 60 days after defect has been detected.

### **How to obtain warranty service?**

For more information on the procedure to follow to obtain warranty service, visit our website at [www.e-swin.com](http://www.e-swin.com).

### **Items not covered by this certificate**

This warranty does not cover damage resulting from accident, non-compliant use (including use other than that described in the user manual of the Device), improper use, negligence, abnormal or unreasonable use, including without limitation, commercial or professional use, adding components not supplied with the Device, lost parts or the use of non-approved spare parts or exposing the device to an electrical current other than that indicated. The warranty will be voided if the Device has been tampered with by a person not approved by **E-Swin** and **E-Swin** will not be liable for any damage resulting

from repairs performed by a person not approved by **E-Jwin**. The warranty does not include replacement of used spare parts.

### **Disclaimer of warranty**

The express warranty set forth in paragraph "14. Warranty " replaces all other warranties, including implied warranties of merchantability or fitness for a particular purpose. This warranty shall be the sole and exclusive remedy in case of defect of the device.

### **Limitations**

In no event shall **E-Jwin** or its approved distributors be liable for any indirect, special, punitive, exemplary, incidental or consequential damage (including, but not limited to damage due to loss of use or cost of replacing the equipment), regardless of the legal theory upon which the claim is based, even if they have been informed of the possibility of such damage in cases where the remedies set forth in this warranty fail in their essential purpose or if **E-Jwin** or its approved distributors are liable for damages under a legal theory. In no event shall the liability of **E-Jwin** or its approved distributors, exceed the total price of a unit, according to the rates in effect on the date of the cause of the action. The user assumes all risk and liability for loss, damage or injury caused to the user, its property, or to other people and their property resulting from improper use or inability to use the product, unless the damage was caused by a negligence from **E-Jwin**. This warranty gives you specific legal rights, to which other rights that vary from country to country may be added. Some countries do not allow the exclusion or limitation of incidental or consequential damages or implied warranties and consequently the limitations and disclaimer above may not apply to you.

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